


DYSSUCTION® SUCTION DEVICE S

REF 67700



REF 67717 (optional)

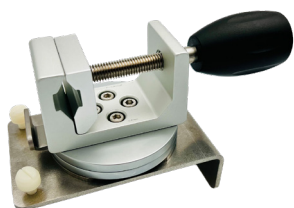

- ▶ Portable suction device for subglottic suction
- ▶ Suction bag with lid for hygienic and germ-reduced replacement
- ▶ Integrated gelling agent for binding secretions
- ▶ Integrated activated charcoal filter for odour reduction in device
- ▶ Digital display in the languages German, English, French, Spanish, Italian, Dutch
- ▶ Individually adjustable
 - Vacuum: -60 mbar to -300 mbar (in 10 mbar increments)
 - Suction time: 10 - 60 sec (setting in 1-second increments)
 - Interval duration: 3 - 60 min (setting in 1-minute increments)

SCOPE OF DELIVERY


DYSSUCTION® SUCTION DEVICE S REF 67700

1 x suction device	
1 x "Bag" external container (1000 ml)	REF 67710
1 x "Bag" external container holder	REF 67712
1 x "OneWay" suction bag (1000 ml)	
1 x disposable secretion tube with stepped connector (non-sterile)	
1 x connecting tube for "OneWay" secretion container system	REF 67713
1 x power supply unit	REF 67912
1 x quick reference guide	
1 x Instructions for Use	

Optional available:

1 x "Bag" external container (1000 ml)	REF 67710
60 x "OneWay" suction bag (1000 ml)	REF 67711
1 x "Bag" external container holder	REF 67712
1 x connecting tube for "OneWay" secretion container system	REF 67713
1 x disposable container system (outer container, suction bag, holder for outer container, connecting tube, disposable suction tube with stepped connector)	REF 67714
10 x disposable secretion tubes with stepped connector (non-sterile)	REF 100712-1
1 x replacement set double filter system DFS®	REF 67716
1 x holder variable DYSSUCTION® Suction Device S	REF 67717

DYSSUCTION® SUCTION DEVICE S

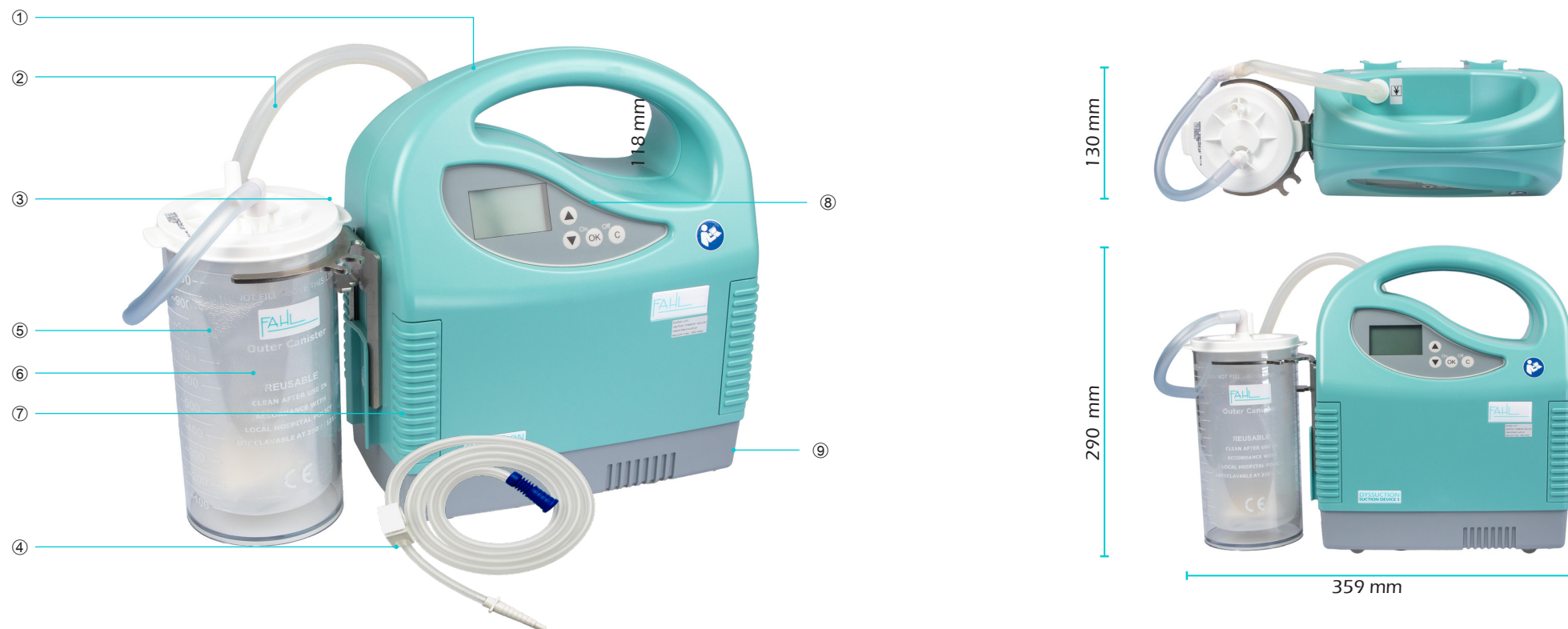

DEVICE OVERVIEW	REF
① Carrying handle	
② Connecting tube for "OneWay" secretion container system	67713
③ Connection suction tube	
④ Disposable secretion tube with stepped connector (non-sterile) (VE 10)	100712-1
⑤ External container "Bag"	67710
⑥ "OneWay" suction bag (1000 ml) (VE 60)	67711
⑦ Holder external container "Bag"	67712
⑧ Control panel	
⑨ Connection for plug-in mains unit	

TECHNICAL DATA	
Vacuum	-60 mbar to -300 mbar (in 10 mbar increments)
Suction time	10 - 60 sec (in 1-second increments)
Interval duration	3 - 60 min (in 1-minute increments)
Unit suction capacity	max. 8 l/min
Dimensions device (H x W x D)	290 x 359 x 130 mm
Weight	2.2 kg
Operating time	Mains operation: permanent operation Battery operation: approx. 18 h depending on the load on the motor
Charging time	6 – 7 h with empty battery
Sound emission	35 dB (A)
Mains voltage mains unit	100 – 240 V AC primary at 50-60 HZ